

# Soluble Sugar Microplate Assay Kit User Manual

Catalog # CAK1186

(Version 1.2A)

Detection and Quantification of Soluble Sugar content in Serum, Plasma, Cell culture, Tissue extracts, Cell lysate, Other biological fluids Samples.

For research use only. Not for diagnostic or therapeutic procedures.



I. INTRODUCTION	2
II. KIT COMPONENTS	3
III. MATERIALS REQUIRED BUT NOT PROVIDED	3
IV. SAMPLE PREPARATION	4
V. ASSAY PROCEDURE	5
VI. CALCULATION	6
VII. TYPICAL DATA	7
VIII. TECHNICAL SUPPORT	7
IX NOTES	7



#### I. INTRODUCTION

Soluble sugar is an important osmotic regulator in plants and is closely related to stress resistance. The stress of plant carbohydrates on low temperature, drought and other stress conditions showed an increase in soluble carbohydrate content of plants. In addition to its important role in metabolism, soluble carbohydrates such as glucose and sucrose help regulate many developmental and physiological processes in plants.

Soluble Sugar Microplate Assay Kit is designed to directly measure soluble sugar content in a variety of samples. In the assay, soluble sugar reacts with anthrone. The product is determined at 620nm, is directly proportional to the soluble sugar concentration in the sample.



#### II. KIT COMPONENTS

Component	Volume	Storage
96-Well Microplate	1 plate	
Dye Reagent	Powder x 1	4 °C
Dye Reagent Diluent	15 ml x 1	4 °C
Standard	Powder x 1	4 °C
Plate Adhesive Strips	3 Strips	
Technical Manual	1 Manual	

#### Note:

Standard: add 1 ml distilled water to dissolve before use, then add 20  $\mu$ l into 980  $\mu$ l distilled water, the concentration will be 200  $\mu$ g/ml.

Dye Reagent: add 15 ml Dye Reagent Diluent to dissolve before use, store at 4 °C.

**Dye Reagent Diluent**: be careful, it is strong acid.

# III. MATERIALS REQUIRED BUT NOT PROVIDED

- 1. Microplate reader to read absorbance at 620 nm
- 2. Distilled water
- 3. Pipettor, multi-channel pipettor
- 4. Pipette tips
- 5. Mortar
- 6. Centrifuge
- 7. Timer



# IV. SAMPLE PREPARATION

# 1. For tissue samples

Weigh out 0.1 g tissue, homogenize with 1 ml distilled water, transfer it into microcentrifuge tubes, incubate at 80 °C for 30 minutes, centrifuged at 10000g for 10 minutes, take the supernatant into a new centrifuge tube for detection.

# 2. For liquid samples

Detect directly.



#### V. ASSAY PROCEDURE

Add following reagents into the microplate:

Reagent	Sample	Standard	Blank
Sample	50 μΙ		
Standard		50 μΙ	
Distilled water			50 μΙ
Dye Reagent	150 μΙ	150 μΙ	150 μΙ

Mix, put the microplate into the oven, incubate at 90 °C for 15 minutes, when cold record absorbance measured at 620 nm.

#### Note:

- 1) Perform 2-fold serial dilutions of the top standards to make the standard curve.
- 2) The concentrations can vary over a wide range depending on the different samples. For unknown samples, we recommend doing a pilot experiment & testing several doses to ensure the readings are within the standard curve range.
- 3) Reagents must be added step by step, can not be mixed and added together.



#### VI. CALCULATION

# 1. According to the volume of sample

SS (
$$\mu$$
g/mI) = (C<sub>Standard</sub> × V<sub>Standard</sub>) × (OD<sub>Sample</sub> - OD<sub>Blank</sub>) / (OD<sub>Standard</sub> - OD<sub>Blank</sub>) / V<sub>Sample</sub>  
= 200 × (OD<sub>Sample</sub> - OD<sub>Blank</sub>) / (OD<sub>Standard</sub> - OD<sub>Blank</sub>)

# 2. According to the weight of sample

SS (
$$\mu$$
g/g) = (C<sub>Standard</sub> × V<sub>Standard</sub>) × (OD<sub>Sample</sub> - OD<sub>Blank</sub>) / (OD<sub>Standard</sub> - OD<sub>Blank</sub>) / (V<sub>Sample</sub> × W/ V<sub>Assay</sub>)
$$= 200 \times (OD_{Sample} - OD_{Blank}) / (OD_{Standard} - OD_{Blank}) / W$$

C<sub>Standard</sub>: the standard concentration, 200 µg/ml;

W: the weight of sample, g;


 $V_{Standard}$ : the volume of standard, 50 µl;

 $V_{Sample}$ : the volume of sample, 50 µl.



# VII. TYPICAL DATA

The standard curve is for demonstration only. A standard curve must be run with each assay.

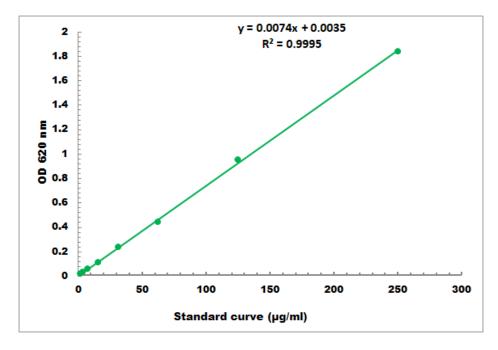

Detection Range:  $2 \mu g/ml - 200 \mu g/ml$ 

# VIII. TECHNICAL SUPPORT

For troubleshooting, information or assistance, please go online to www.cohesionbio.com or contact us at techsupport@cohesionbio.com

# IX. NOTES