

Product Data Sheet

Anti-Ki67 Antibody

Catalog # Source Reactivity Applications

CPA9691 Mouse H IH, IF

Description Mouse monoclonal antibody to Ki67

Immunogen KLH-conjugated synthetic peptide encompassing a sequence within human Ki67. The

exact sequence is proprietary.

Purification The antibody was purified by immunogen affinity chromatography.

Specificity Recognizes endogenous levels of Ki67 protein.

Clonality Monoclonal

Conjugation

Form Mouse IgG1. Liquid in PBS containing 50% glycerol, 0.2% BSA and 0.01% sodium

azide.

Dilution IH (1/100 - 1/300), IF (1/100 - 1/500)

Gene Symbol MKI67

Alternative Names Antigen KI-67

Entrez Gene 4288 (Human)

SwissProt P46013 (Human)

Storage/Stability Shipped at 4°C. Upon delivery aliquot and store at -20°C for one year. Avoid

freeze/thaw cycles.

Application key: E- ELISA, WB- Western blot, IH- Immunohistochemistry, IF- Immunofluorescence, FC- Flow cytometry, IC- Immunocytochemistry, IP- Immunoprecipitation, ChIP- Chromatin Immunoprecipitation, EMSA- Electrophoretic Mobility Shift Assay, BL- Blocking, SE- Sandwich ELISA, CBE- Cell-based ELISA, RNAi- RNA interference Species reactivity key: H- Human, M- Mouse, R- Rat, B- Bovine, C- Chicken, D- Dog, G- Goat, Mk- Monkey, P- Pig, Rb-Rabbit, S- Sheep, Z- Zebrafish

COHESION BIOSCIENCES LIMITED

WEB ORDER SUPPORT CUSTOM
www.cohesionbio.com order@cohesionbio.com techsupport@cohesionbio.com custom@cohesionbio.com










Immunohistochemical analysis of Ki67 staining in human classical hodgkin's lymphoma formalin fixed paraffin embedded tissue section. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.0). The section was the

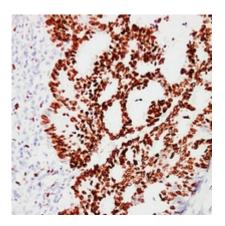

Immunohistochemical analysis of Ki67 staining in human colon carcinoma formalin fixed paraffin embedded tissue section. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.0). The section was then incubated with the antibody at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

Application key: E- ELISA, WB- Western blot, IH- Immunohistochemistry, IF- Immunofluorescence, FC- Flow cytometry, IC-Immunocytochemistry, IP- Immunoprecipitation, ChIP- Chromatin Immunoprecipitation, EMSA- Electrophoretic Mobility Shift Assay, BL- Blocking, SE- Sandwich ELISA, CBE- Cell-based ELISA, RNAi- RNA interference Species reactivity key: H- Human, M- Mouse, R- Rat, B- Bovine, C- Chicken, D- Dog, G- Goat, Mk- Monkey, P- Pig, Rb-Rabbit, S- Sheep, Z- Zebrafish